



health
care
professional

Prizma GmbH
GERMANY
40213 Düsseldorf, Graf-Adolf-Platz, 15
Tel. +49 211 97632861
e-mail: s.katanic@prizma.gmbh
www.prizmade.de

Prizma Ltd
SERBIA
34000 Kragujevac, N° 8, Kumanovska st.
Tel: +381 34 330 200
Fax: +381 34 336 146, 336 147
e-mail: sales@prizma.rs
www.prizma.rs





OUR STORY

1990
year that we were founded

Production Experience
25 years

2019, foundation of Prizma GmbH,
based in Germany

EU, USA, Indonesia, Australia
our export markets

From 1990 Omron Healthcare Distributor
Serbia, Bosnia & Herzegovina, Montenegro
and North Macedonia markets


Cooperation with
Rwth Aachen University, Serbian Academy
Of Science

CE  CE 1304
our products certificates

ISO 9001:2015, EN ISO 13485:2016, ISO 14001:2015
our management system standards



breathe care

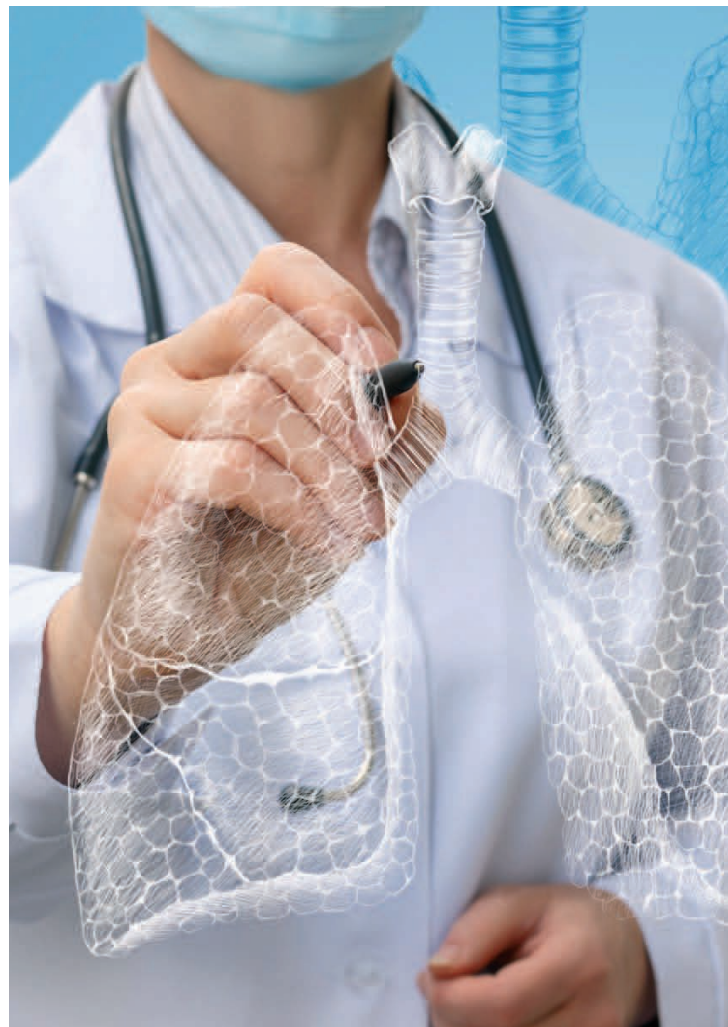



MEDICAL DEVICE

EFFICIENT OXYGEN DELIVERY SYSTEM

PRIZMA O2, portable oxygen concentrator, provides a high oxygen concentration to patients, with oxygen needs.

The device is delivering medical-grade oxygen through the nasal cannula, while removing nitrogen.

When a patient needs oxygen, a system for effective and efficient oxygen delivery is essential.

PRIZMA O2 is a lightweight device that comes with a battery and can be used at home, at work or while travelling.

PRIZMA O2 delivers oxygen through a pulse dose using intelligent technology.

PRIZMA O2 has simple control functions and an LCD screen.

PRIZMA O2

CE 1304

- **Oxygen delivery**
on the go
- **Lightweight**
and quiet
- **Easy to use**
and easy to maintain



Technical Data

- Oxygen Concentration: 90-96%
- Flow type: pulse
- Levels of flow settings: 6
- Flow rate: 210-1260 ml/min
- Max limited pressure: 1.5 bar \pm 25%
- Battery life: up to 6 h, at 210 ml/min
- Battery charging time: 2.5 h
- Noise level: \leq 60 dB
- Power Supply: 100-240 V AC, 50-60 Hz; 24 V DC, 145 VA
- Warm up time: 2 min
- Dimensions: 187 \times 88 \times 217 mm
- Weight: 2.16 kg



DEVELOPED & MADE IN SERBIA

CONTENTS OF THE KIT:

PRIZMA O2 oxygen concentrator, 8-core battery, AC/DC adapter, nasal cannula, carrying bag and instruction manual



MEDICAL DEVICE

PROFESSIONAL ULTRASONIC NEBULIZERS

Prizma is producing the own line of professional ultrasonic nebulizers for medical facilities; devices are simple to use and easy to maintain. Our professional ultrasonic nebulizers are ideal for continuous multipatient nebulization.

On the device's touch screen it is possible to make settings such as: air volume, nebulization rate, start time and duration of inhalation therapy. Screen is readable even in the dark, thanks to the backlight.

Our product range consists of:

- PROFISONIC** – basic professional ultrasonic nebulizer
- PROFISONIC H** – professional ultrasonic nebulizer with integrated aerosol heating system



DEVELOPED & MADE IN SERBIA

PROFI SONIC

CE 1304

- **Large number of patients**
continuous inhalations
- **Touch screen**
for parameters setting
- **In every position**
nebulization



Technical Data

- Nebulization rate: 0.0-3.0 ml/min
- Air volume: max 20 l/min
- Particle size: 0.5-5.0 microns
- Capacity of medication cup: 5-150 ml
- Low drug residue: 1 ml
- Capacity of water tank: 300 ml
- 100% silicone hose:
ultra-thin silicone hose can be sterilized in an autoclave
- Ultrasonic frequency: 1.7 MHz
- Noise level: <35 dB (A)
- Operating voltage/frequency: 100-240 V AC, 50/60 Hz
- Power consumption: 55 W
- Dimensions: 260 × 200 × 250 mm
- Weight: 1.8 kg

CONTENTS OF THE KIT:

PROFI SONIC nebulizer, silicone inhalation hose (90 cm), mouthpiece, 3 medication cups (one already installed), adult mask, child mask, bacterial filter, power cord and instruction manual.

PROFI SONIC H

CE 1304

- **Heated nebulization**
continuous inhalations
- **Adjusting the heat level**
of nebulization
- **Touch screen**
for parameters setting

Technical Data

- Nebulization rate: 0.0-3.0 ml/min
- Air volume: max 20 l/min
- Particle size: 0.5-5.0 microns
- Capacity of medication cup:
5-150 ml
- Low drug residue: 1 ml
- Capacity of water tank: 300 ml
- 100% silicone hose:
ultra-thin silicone hose can be sterilized in an autoclave
- Ultrasonic frequency: 1.7 MHz
- Noise level: <35 dB (A)
- Operating voltage/frequency: 100-240 V AC, 50/60 Hz
- Power consumption: 90 VA
- Dimensions: 260 × 200 × 250 mm
- Weight: 1.8 kg

CONTENTS OF THE KIT:

PROFI SONIC H nebulizer, silicone inhalation hose (90 cm),
3 medication cups (one already installed), bacterial filter, power cord
and instruction manual.

Optional accessories are available.

OPTIONAL ACCESSORIES

PROFI SONIC / PROFI SONIC H

Table stand, adjustable height
Made of stainless steel
Catalogue code: PSDS



PROFI SONIC / PROFI SONIC H


Movable stand, adjustable height
Made of stainless steel
Catalogue code: PSTS



PROFI SONIC / PROFI SONIC H
Medication supply set
Catalogue Code: PSMS
Material: PP (polypropylene)
Single pack quantity: 1 piece



PROFI SONIC / PROFI SONIC H
Inhalation hose
Catalogue Code: PSIH
Material: Silicone
Single pack quantity: 1 piece

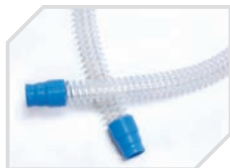


PROFI SONIC / PROFI SONIC H
Bacterial filter
Catalogue Code: PSBF
Material: Polypropylene
Single pack quantity: 1 piece



PROFI SONIC / PROFI SONIC H
Medication cup
Catalogue Code: PSMC
Material: PET
Single pack quantity: 5 pieces



PROFI SONIC
Mouthpiece
Catalogue Code: PSMP
Material: PC
Single pack quantity: 5 pieces

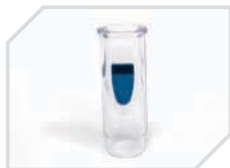

PROFI SONIC
Child mask
Catalogue Code: PSCM
Material: Silicone
Single pack quantity: 1 piece



PROFI SONIC
Adult mask
Catalogue Code: PSAM
Material: Silicone
Single pack quantity: 1 piece





MEDICAL DEVICE

COMPRESSOR NEBULIZER

Prizma is producing powerful compressor nebulizer with unique nebulizer kit that allows setting of nebulization rate of 0.3-0.45 ml/min, providing efficient drug delivery to the smallest parts of the airways and making the therapeutic effect of the administered drug even greater.

During nebulization the users are protected of inhaling dust and other impurities from the air, thanks to the filter that is supplied in our compressor nebulizer.

Our compressor nebulizer is intended for use in the home environment and medical facilities, and patients could be adults and kids.



DEVELOPED & MADE IN SERBIA

PRIZJet

CE 1304



- **Setting of aerosol rate with unique nebulizer kit**
- **Strong compressor**
- **For home & medical facilities use**

Technical Data

- Particle size: <5 microns
- Nebulization rate: 0.3-0.45 ml/min
- Capacity of medication cup: 2-8 ml
- Low drug residue: 0.64 ml
- Compressor maximum pressure: 2.5 bar
- Compressor maximum air flow: 8 l/min
- Noise level: <60 dB (A)
- Operating voltage/frequency: 230 V/50 Hz
- Power consumption: 120 VA
- Dimensions: 180 × 140 × 110 mm
- Weight: 1.3 kg (compressor only)

CONTENTS OF THE KIT:

PRIZJet nebulizer, the power cord, inhalation hose (length 150 cm), nebulizer kit (medication cup), adult mask, child mask, mouthpiece, air filters (5 pcs + 1 assembled), storage and carrying case pouch and instruction manual.

Optional accessories available for different types of inhalation therapy.



wellness to health





MEDICAL DEVICE

DRY SALT AEROSOL GENERATORS

Prizma is producing the own line of dry salt aerosol generators with EC and UL certificate for medical facilities and wellness / SPA centers.

On the device's touch screen or via the Android / Windows / iOS smart device is possible to make settings of device, such as salt consumption, start time, duration of halotherapy, dosage profile (device is offering 10 dosage profiles). Screen is readable even in the dark, thanks to the backlight.

Materials, which are in contact with salt, are made of stainless steel and plastic, that ensures long life of Prizma dry salt aerosol generators.

Our product range consists of:

- multiSALT+** – wellness / SPA device with UL/EC certificate for small salt rooms, salt cabins and tent
- PRIZSALT+** – medical device with EC certificate
- PRIZSALT+S** – wellness / SPA device with UL/EC certificate
- homeSALT** – home device with UL/EC certificate to be used also in salt room, tent, cabin or in outdoor spaces
- nanoSALT** – home device with UL/EC certificate to be used also in salt room, tent, cabin or in outdoor spaces



DEVELOPED & MADE IN SERBIA

multiSALT+



- **Wellness/SPA device**
- **Wireless management**
Prizma Link app
- **Up to 9 programs**
personalization



Technical Data

- Touch screen operated
- Personalization of up to 9 programs: duration of inhalation and salt dosing
- Capacity of the salt container: max 35 g
- Salt consumption: 0-100 g/h (adjustable)
- Operating voltage/frequency: 100-240 V AC, 50/60 Hz
- Power consumption: 18 W
- Dimensions: 225 × 218 × 395 mm
- Weight: 8 kg

CONTENTS OF THE KIT:

multiSALT+ dry salt aerosol generator, cleaning equipment, power cord, deflector, extension tube, software and instruction manual.



PRIZSALT+

CE 1304

- **Medical device**
- **Wireless management**
Prizma Link app
- **Up to 9 programs**
personalization



PRIZSALT+S

UL CE

- **Wellness/SPA device**
- **Wireless management**
Prizma Link app
- **Up to 9 programs**
personalization



Technical Data

- Touch screen operated
- Personalization of up to 9 programs:
duration of inhalation and Salt dosing
- Capacity of salt dispenser: max 150 g
- Salt consumption: 0-500 g/h (adjustable)
- Operating voltage/frequency:
100-240 V AC, 50/60 Hz
- Power Consumption: 35 W
- Dimensions: 270 × 260 × 480 mm
- Weight: 13.5 kg

CONTENTS OF THE KIT:

PRIZSALT+ dry salt aerosol generator, cleaning equipment, power cord, deflector, extension tube, software and instruction manual.



Technical Data

- Touch screen operated
- Personalization of up to 9 programs:
duration of inhalation and Salt dosing
- Capacity of the salt dispenser: max 150 g
- Salt consumption: 0-500 g/h (adjustable)
- Operating voltage/frequency:
100-240 V AC, 50/60 Hz
- Power consumption: 35 W
- Dimensions: 270 × 260 × 480 mm
- Weight: 13.5 kg

CONTENTS OF THE KIT:

PRIZSALT+S dry salt aerosol generator, cleaning equipment, power cord, deflector, extension tube, software and instruction manual.



homeSALT



- **New approach in halotherapy**
- **Home device for indoor and outdoor spaces**
- **Remote control**
4 programs of operation

Technical Data

- Remote control: Operate from a distance
- Programs of operation: 4 programs of operation
- Led indicator: Light indication of device status
- Operating voltage: 24 V
- AC/DC adapter MPU60B-108 :
100-240 V AC, 47-63 Hz, 1.6-0.7 A
24 V DC, 2.6 A max
- Power consumption: 15 W
- El.protection class: Class III
- Salt tank capacity: 20 g maximum
- Dimensions: 113 × 236 × 70 mm
- Net weight: 1 kg

CONTENTS OF THE KIT:

homeSALT dry salt aerosol generator, AC/DC adapter, remote control key fob, brush and instruction manual.



nanoSALT



- **New approach in halotherapy**
- **Small dimension battery operated**
- **Portable**

Technical Data

- 10 inhalations with a fully charged battery
- Battery charge indicator
- Empty battery indicator
- Indicator of blocked motor
- Open cover indicator
- Operating voltage: 6V DC
- AC/DC adapter:
100-240 V AC, 47-63 Hz, 1.6-0.7 A
5 V 3 A DC
- Power consumption: 5 W
- El.protection class: Class III
- Salt tank capacity: 5 g
- Dimensions: 170 × 50 × 60 mm
- Net weight: 0.5 kg

CONTENTS OF THE KIT:

nanoSALT dry salt aerosol generator, AC/DC adapter, brush and instruction manual.





high precision





MEDICAL DEVICE

ANEROID SPHYGMOMANOMETER

Prizma is producing aneroid sphygmomanometer with EC certificate.

100% German mechanism made by HB-Feinmechanik GmbH & Co. KG is the heart of our aneroid sphygmomanometer, that provides accurate measurements and reliable readings.

Its one-tube system guarantees easiest operation while unique designed 180° rotary anti-slip bulb enables both left- and right-hand users to easily use the device.

Bladderless cuffs are biocompatible according to the high criteria of EN ISO 10993 series of standards and produced with high quality TPU coated nylon fabric that are latex-free, not toxic, non phthalate.

It is intended for use in the medical facilities and in home environment.

PRIZMA PA1

CE 1304



- **100% German mechanism** by HB-Feinmechanik
- **180° rotary anti-slip bulb** for left-handed and right-handed users
- **Bladderless cuff** biocompatible (EN ISO 10993)

Technical Data

- Manometer scale range: 0-300 mmHg
- Accuracy: ± 3 mmHg
- Scale diameter: $\varnothing 66$ mm
- Tube system: one-tube system
- Cuff type: bladderless for arm circumference 27.5-36.5 cm
- Cuff material: TPU coated nylon
- Measuring mechanism, German-made by HB-Feinmechanik GmbH & Co. KG

CONTENTS OF THE KIT:

PRIZMA PA1 aneroid sphygmomanometer, bladderless cuff 27.5-36.5 cm, storage and carrying case pouch and instruction manual. Optional cuffs for arm circumferences of 9.0-14.8 cm, 13.8-21.5 cm, 20.5-28.5 cm, 27.5-36.5 cm and 35.5-46.0 cm are available on demand



DEVELOPED & MADE IN SERBIA

long tradition





INDUSTRIAL DEVICE

R&D SCALE ULTRASONIC SPRAY PYROLYSIS SYSTEM

PRIZnano is a large capacity ultrasonic atomizer of liquids. Usage of 3 piezo-electric crystals guarantee a large atomization capacity up to 1.2 l/h, which provides production of large quantity of aerosol that is ideal for independent use (humidification of large spaces) or in combination with other system for the production of nano particles.

PRIZMA ESP nano particle collector use high-voltage electrostatic field to collect ultra-fine particles. The minimum collectable particle is 0.1 nm. Collection efficiency is more than 95%.

It is intended for use in industry, science and research & development centers (in the process of obtaining nano powders), in medicine and food industry.



DEVELOPED & MADE IN SERBIA

PRIZnano

CE



- **Pyrolysis atomizer**
our specialization
- **Resistant to chemicals**
the chamber
- **RWTH Aachen university**
developed with


Technical Data

- 3 piezoelectric crystals, which are glass coated due to resistance to chemicals
- Capacity of atomization: 0-1.2 l/h
- The most common particle size:
 - 5.0 microns - 1.7 MHz
 - 3.0 microns - 2.4 MHz
- Ultrasonic frequency: 1.7 or 2.4 MHz
- Power supply / frequency:
100-240 V AC, 50/60 Hz
- Power consumption: 150 W
- Dimensions: 250 × 250 × 350 mm
- Weight: 6.5 kg

CONTENTS OF THE KIT:

PRIZnano atomizer and instruction manual.

PRIZMA ESP

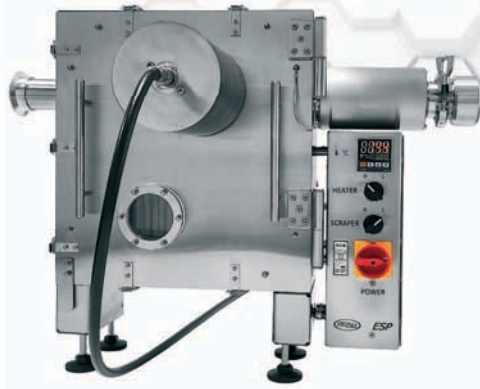

notes

- Easy collection of nano powder
- Easy cleaning



Collected silver nano particles

Technical Data

- Power supply / frequency:
100-240 V AC, 50/60 Hz
- Power consumption: 500 W
- Dimensions: 600 × 550 × 560
- Weight: 32.5 kg

CONTENTS OF THE KIT:

PRIZMA ESP filter and instruction manual.

